




# NewCare Junior Bed

## QUICK REFERENCE SHEET

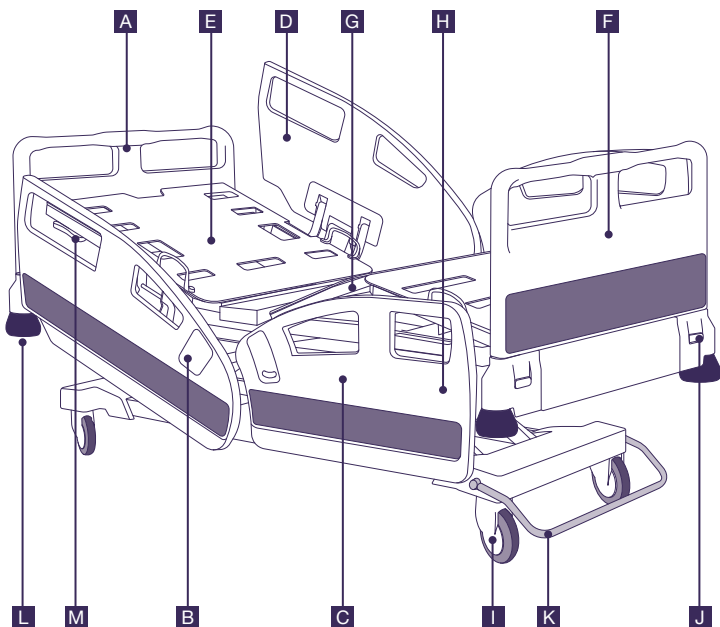

### BED'S COMPONENTS

A	Removable curved headboard	H	Legrest side rails
B	Backrest tilting indicator	I	Double wheels Ø15 cm
C	Trend / Reverse Trend tilting indicator	J	Headboard and footboard unlocking
D	Backrest side rails	K	Braking pedal
E	HPL mattress platform	L	Protection discs
F	Removable curved footboard	M	Emergency CPR
G	Sheet steel central section		

### OTHER FUNCTIONS

#### Side Rails (1):


- Raising the rail: Pull the rail upwards until it is locked in the highest position.
- Lowering the rail: Pull the handle found underneath the rail outwards. The rail will automatically lower.

#### Emergency CPR (2):


- Activate either of the two levers located under the headrest module, accompanying it in turn as it lowers.

#### Centralised braking system (3):



- Braking the bed: Pedal in low position.
- Moving the bed: Pedal in the middle position.
- Activating the 5th or steering wheel: Pedal in raised position.

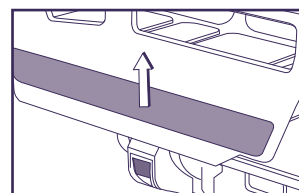


Rising the rails

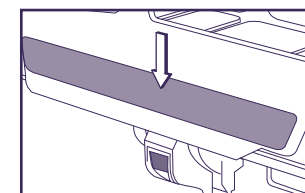


Lowering the rails

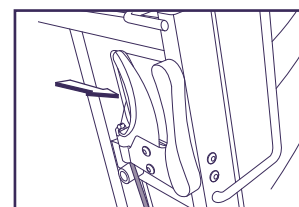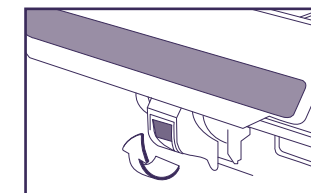


Emergency CPR

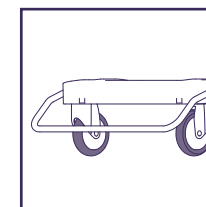

Brake

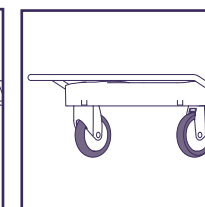

Free movement

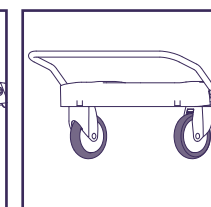

Steering or 5th wheel

### CLEANING

Daily cleaning and disinfection with up to a 1000ppm concentration of chlorine solutions:

- Rails, headboard, footboard.

Cleaning and disinfection after a patient leaving or before a transfer:

- Mattress platform area, chassis of the bed.

Deep cleaning and disinfection every two months:

- Support and lifting area of the bed.

### MAINTENANCE

Always check the correct state of:

- Electrical wiring.
- Wheels.
- Braking mechanisms.

www.pardo.es


TECHNICAL SERVICE CENTRE

(+34) 976 300 033


For more information, please consult the use and maintenance manual.




# CONTROLS

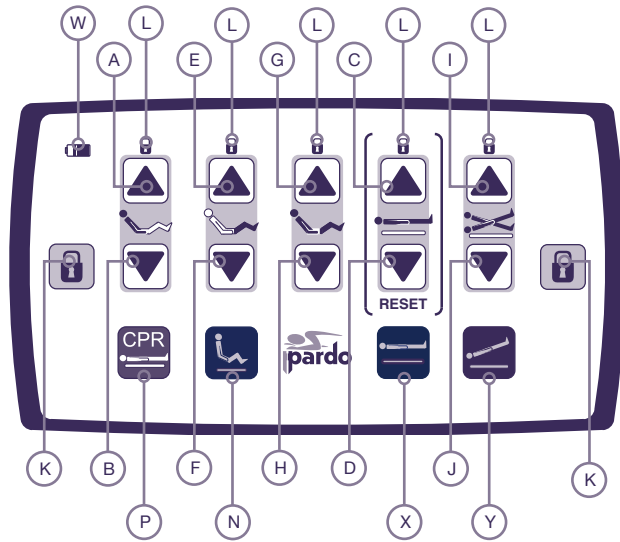


Wired nursing control

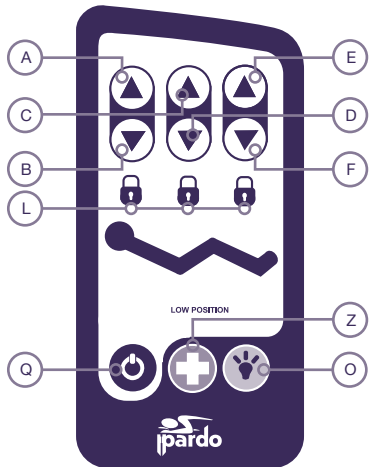


Membrane control w/o locks

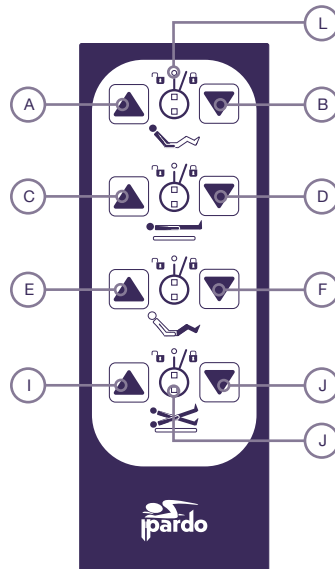

Patient control with locks

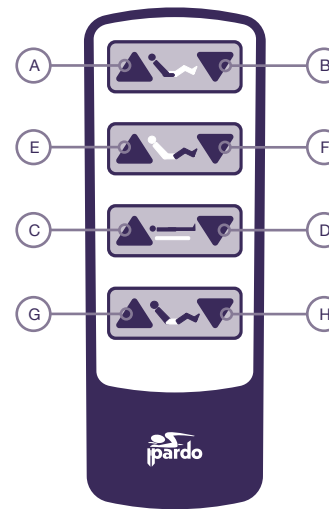

Patient control

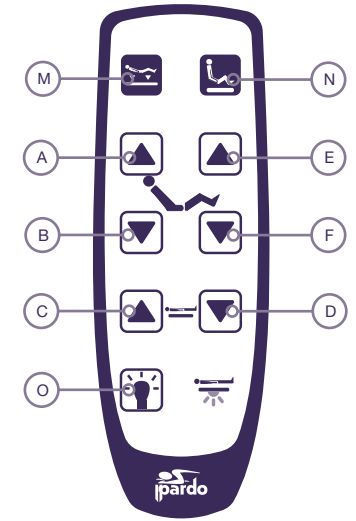

Patient control with light

## INSTRUCTIONS OF USE

A	Rise backrest	H	Lower both sections simultaneously ( <i>autocontour</i> )	O	Light under bed
B	Lower backrest	I	Reverse Trend	P	Emergency CPR
C	Rise mattress platform	J	Trend	Q	Standby button
D	Lower mattress platform	K	Locking movements	W	Battery indicator
E	Rise legrest	L	Indicator Led of locked movements	X	Exploration Position (maximum height)
F	Lower legrest	M	Flat Position (lower all sections simultaneously)	Y	Emergency Trend position
G	Rise both sections simultaneously ( <i>autocontour</i> )	N	Cardiac chair position	Z	Nursing call

Industrias Hidráulicas Pardo S.L. company has the following certifications:



### SAFETY WARNINGS

- The maximum duration of a continuous manoeuvre should not exceed 2 minutes, as damage may be caused to the equipment.
- Before raising and lowering the moving parts of the mattress platform, with the other parts of the bed or accessories, it is essential to check that there are no obstacles between them.
- The Trendelenburg and reverse Trendelenburg movement should always be performed under medical supervision.
- When the bed is in extended position, do not sit or lean on the foot part of the bed, since it may cause a bed dump.
- Do not use CPR bilateral actuating levers to raise the backrest of the bed, it may cause raise motor breakdown and makes a bed dysfunction.

